


# Nano & UV Artificial Teeth Cleaner

Clean your artificial teeth with new technology-Nano-silver & UV light.

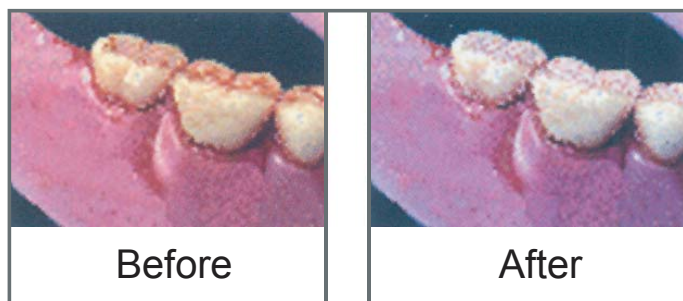

## DOES WASHING ARTIFICIAL TEETH ALWAYS BOTHER YOU?

With the brand-new "Nano & UV Artificial Teeth Cleanser", you can easily clean your smudgy denture. It efficiently clean your artificial teeth or jewelry with nano-silver and latest Japan technology - UV light sterilization.

### Product Features:

- \* Through the latest UV light sterilization and nano-silver bacteriostasis to reach 100% spotless effect without residual chemical.
- \* Clean dentures with the high power of 8,000 times per minute vibration.
- \* Build-in power save timer to shut the power automatically after 5 minutes when you turn it on.
- \* Small size is perfect for home-use (by AC adapter) and travel-use (Battery)
- \* UV light sterilization and nano-silver mixture container can deodorize, sterilize and wash at one step.

### Product Specification:

- \* Dimension : 120 (L) x 130 (W) x 105(D) (mm) / 4.72 (L) x 5.12 (W) x 4.13 (D) (in)
- \* Weight : 370g
- \* Vibration : 8,000 times per minute
- \* Power Supply : Linear AC adapter 110V~240V.

### Illustration:

1. Install batteries in the unit or plug the adapter in.
2. Take out the container and pour onto it reach 80% full.
3. Put the container back properly & tightly on the cleaner.
4. Set down the denture or jewelry and then choose either one function.  
(Washing & UV or UV only). It powers off automatically after 5 minutes when you turn it on.
5. In case you want to stop working in advance, press the function button for 3 seconds.
6. When the washing is done, take it out and pour out water. ( The unit is waterproof. )



DIGIMAX INNOVATIVE PRODUCTS LTD.

2F, No.196, Sec. 2, Zhong-Xing Road,  
Hsin-Tien City, Taipei Hsien 231, Taiwan R.O.C.  
886-2-8665-0966 tel ext:866 886-2-8665-0957 fax  
E-mail: monika@digimaxproducts.com  
<http://www.digimaxproducts.com>